


MACHEREY-NAGEL

# User manual

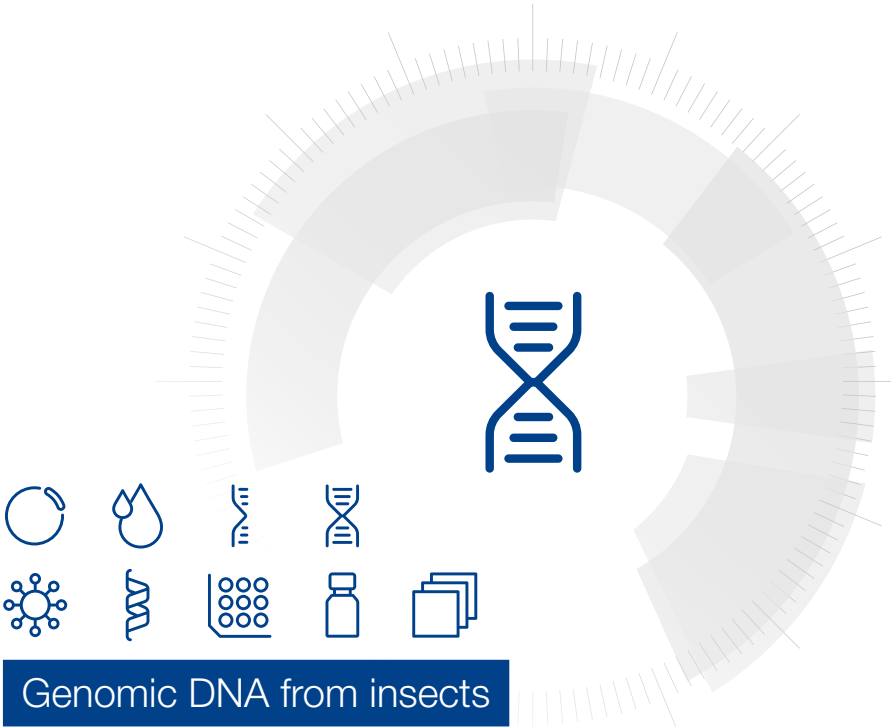


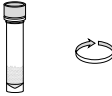
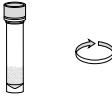
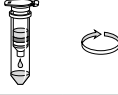
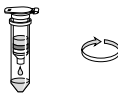
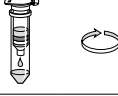
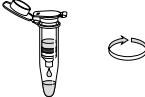
■ NucleoSpin® DNA Insect

October 2023 / Rev. 04

# Genomic DNA from insects

## Protocol at a glance (Rev. 03)

### NucleoSpin® DNA Insect

1 Prepare sample		Place < 40 mg insect material (wet weight) in NucleoSpin® Bead Tube Type D  100 µL Elution Buffer BE						
2 Lyse sample		40 µL Buffer MG  10 µL Liquid Proteinase K  Agitate on a swing mill, MN BeadTube Holder or similar device 0.5 – 15 min  11,000 × g, 30 s						
3 Adjust binding conditions		600 µL Buffer MG  11,000 × g, 30 s						
4 Bind DNA		Load 500 – 600 µL sample on NucleoSpin® DNA Insect Column  11,000 × g, 30 s						
5 Wash silica membrane		<table border="0"> <tbody> <tr> <td style="background-color: black; color: white; padding: 2px;"><b>1<sup>st</sup></b></td> <td>500 µL BW</td> <td>11,000 × g, 30 s</td> </tr> <tr> <td style="background-color: black; color: white; padding: 2px;"><b>2<sup>nd</sup></b></td> <td>500 µL B5</td> <td>11,000 × g, 30 s</td> </tr> </tbody> </table>	<b>1<sup>st</sup></b>	500 µL BW	11,000 × g, 30 s	<b>2<sup>nd</sup></b>	500 µL B5	11,000 × g, 30 s
<b>1<sup>st</sup></b>	500 µL BW	11,000 × g, 30 s						
<b>2<sup>nd</sup></b>	500 µL B5	11,000 × g, 30 s						
6 Dry silica membrane		11,000 × g, 30 s						
7 Elute DNA		100 µL Elution Buffer BE  RT, 1 min  11,000 × g, 30 s						

---

## Table of contents

1	Components	4
1.1	Kit contents	4
1.2	Reagents, consumables, and equipment to be supplied by user	5
1.3	About this user manual	5
2	Product description	6
2.1	The basic principle	6
2.2	Kit specifications	6
2.3	Handling, preparation, and storage of starting materials	6
2.4	Lysis and disruption of sample material	7
2.4.1	Disruption device	7
2.4.2	Type of bead tube	7
2.4.3	Time and frequency of disruption	8
2.5	Elution procedures	9
3	Storage conditions and preparation of working solutions	10
4	Safety instructions	11
4.1	Disposal	11
5	Protocols	12
5.1	Protocol for fresh, frozen, dried, and ethanol preserved insect samples	12
5.2	Protocol for purification of DNA from hard to lyse bacteria in insect samples	14
6	Appendix	15
6.1	Troubleshooting	15
6.2	Ordering information	17
6.3	Product use restriction/warranty	18

# 1 Components

## 1.1 Kit contents

NucleoSpin® DNA Insect		
REF	10 preps 740470.10	50 preps 740470.50
Lysis Buffer MG	10 mL	38 mL
Wash Buffer BW	6 mL	30 mL
Wash Buffer B5 (Concentrate)*	6 mL	6 mL
Elution Buffer BE**	13 mL	30 mL
Liquid Proteinase K	120 µL	600 µL
NucleoSpin® Bead Tubes Type D	10	50
NucleoSpin® DNA Insect Columns (light green rings)	10	50
Collection Tubes (2 mL)	20	100
User manual	1	1

\* For preparation of working solutions and storage, see section 3.

\*\*Composition of Elution Buffer BE: 5 mM Tris/HCl, pH 8.5

## 1.2 Reagents, consumables, and equipment to be supplied by user

### Reagents

- 96–100 % ethanol (for preparation of Wash Buffer B5)

### Consumables

- 1.5 mL or 2 mL microcentrifuge tubes for sample preparation and elution
- Disposable tips

### Equipment

- Manual pipettors
- Centrifuge for microcentrifuge tubes
- Vortex mixer
- Personal protection equipment (lab coat, gloves, goggles)
- Sample disruption device:  
The MN Bead Tube Holder (REF 740469, see ordering information, section 6.2) is recommended to be used in combination with the Vortex-Genie® 2 for cost efficient and convenient disruption of lipid tissue samples. The Vortex Adapter (MoBio) for Vortex-Genie® 2 X is also suitable.

Alternatively, a swing mill can be used considering precautions of section 2.4.3 (e.g., mixer mill MM200, MM300, MM400 (Retsch®)).

**WARNING:** *The use of other disruption devices like FastPrep® System (MP-Biomedicals), Precellys® (Bertin Technologies), MagNA™ Lyser (Roche), TissueLyser (QIAGEN), Bullet Blender® (Next Advance), Mini-Beadbeater™ (Biospec Products), Speed Mill (Analytik Jena), or similar devices might cause bead tube destruction. Such disruption devices can cause high mechanical stress on the bead tubes. Depending on bead tube type and content (beads like steel beads, liquid volume, sample type), especially high frequency of shaking and / or long shaking duration can cause destruction of the bead tubes. **If using such a disruption device, it is the responsibility of the user to perform initial stability tests to ensure stability of bead tubes during the individual experimental setup (e.g., intensity of agitation).** See also section 2.4.3!*

## 1.3 About this user manual

It is strongly recommended for first time users to read the detailed protocol sections of the **NucleoSpin® DNA Insect** kit before using this product. Experienced users, however, may refer to the Protocol at a glance instead. The Protocol at a glance is designed to be used only as a supplemental tool for quick referencing while performing the purification procedure.

All technical literature is available online at [www.mn-net.com](http://www.mn-net.com).

Please contact Technical Service regarding information about any changes to the current user manual compared with previous revisions.

## 2 Product description

### 2.1 The basic principle

The **NucleoSpin® DNA Insect** kit is designed for efficient isolation of genomic DNA from insect samples. DNA can be isolated from a wide variety of samples, such as fresh, frozen, dried and ethanol preserved insects, e.g., fruit fly (*Drosophila melanogaster*), house cricket (*Acheta domesticus*), field cricket (*Gryllus assimilis*), mealworm (*Tenebrio molitor*), non-biting midge larvae (*Chironomidae*), and mosquito (*Culicidae*).

Insects can be difficult to lyse due to their strong, chitin reinforced cell walls.

The **NucleoSpin® DNA Insect** kit combines enzymatic lysis by utilizing mechanical disruption of cell walls with the NucleoSpin® Bead Tubes. The NucleoSpin® Bead Tubes can be used in combination with the MN Bead Tube Holder (REF 740469) and the Vortex-Genie® 2. They are also compatible with other disruptive devices (see section 2.4.1). High DNA yields can be obtained with the NucleoSpin® Bead Tubes from a large variety of sample types – enabling the procedure to be convenient, fast, and easy. Alternative bead tube types can be ordered separately for selected sample types (see section 2.4.2 for recommendations).

### 2.2 Kit specifications

#### Kit specifications at a glance

Parameter	NucleoSpin® DNA Insect
Technology	Silica membrane technology
Format	Mini spin column
Sample material	Fresh, frozen, dried, and ethanol preserved insect samples
Sample amount	Up to approx. 40 mg wet weight
Typical yield	Varies by sample and disruption device. Up to 25 µg DNA can be obtained.
$A_{260} / A_{280}$	1.7 – 1.9
Elution volume	25 – 200 µL
Preparation time	35 min/6 preps
Binding capacity	60 µg
Use	For research use only

### 2.3 Handling, preparation, and storage of starting materials

Fresh, frozen, dried, and ethanol preserved insect samples can be used. Make sure not to use more than 40 mg starting material.

## 2.4 Lysis and disruption of sample material

In order to obtain optimal DNA yields, a complete disruption of the sample material is essential. The efficiency of sample disruption depends on the parameters listed below and suggestions for optimization are outlined in the subsequent sections.

### 2.4.1 Disruption device

The following devices are compatible with NucleoSpin® Bead Tubes. Please check whether NucleoSpin® Bead Tubes can be accommodated by the available tube adapters prior to starting the procedure.

- MN Bead Tube Holder in combination with the Vortex-Genie® 2 (recommended).
- Mixer mill MM200, MM300, MM400 (Retsch®) (suitable).

*If other disruption devices (section 1.2) are intended to be used, consider section 2.4.2 and WARNING note in section 2.4.3!*

### 2.4.2 Type of bead tube

Bead type, disruption time, and frequency / speed must be optimized for a given sample to obtain maximal DNA yield and quality.

#### Type of bead tube

- NucleoSpin® Bead Tubes Type D (3 mm steel beads; included in NucleoSpin® DNA Insect kits)  
Recommended for insect samples

For more information see here:

**<https://www.mn-net.com/de/BeadTubeOverview-BeadTubeD>**

Other types of bead tubes are available for other applications:

- NucleoSpin® Bead Tubes Type A (0.6–0.8 mm ceramic beads)  
Recommended for soil and sediment (included in NucleoSpin® Soil, see ordering information, section 6.2).
- NucleoSpin® Bead Tubes Type B (40–400 µm glass beads)  
Recommended for gram positive and -negative bacteria (included in NucleoSpin® Microbial DNA, see ordering information, section 6.2).
- NucleoSpin® Bead Tubes Type C (1–3 mm corundum)  
Recommended for yeast (see ordering information, section 6.2).
- NucleoSpin® Bead Tubes Type E (combination of 3 mm steel beads and 40–400 µm glass beads)  
Recommended for difficult to lyse tissue containing gram positive bacteria (see ordering information, section 6.2).
- NucleoSpin® Bead Tubes Type F (1–3 mm corundum and 3 mm steel beads)  
Use only with MN Bead Tube Holder! (see ordering information, section 6.2)

## 2.4.3 Time and frequency of disruption

The following recommendations have been established for the MN Bead Tube Holder in combination with a Vortex-Genie® 2 or a Retsch® Swingmill MM300 operating at highest frequency (30 Hertz). For using other disruption devices, and other sample materials, time and frequency have to be optimized.

Sample material	Disruption device	Disruption time	Speed / intensity / frequency
Fresh, frozen, dried, and ethanol preserved insects, e.g., <i>Drosophila melanogaster</i>	MN Bead Tube Holder in conjunction with a Vortex-Genie® 2	approx. 20 min	full speed
Fresh, frozen, dried, and ethanol preserved insects, e.g., <i>Drosophila melanogaster</i>	Mixer mill (Retsch®)	approx. 0.5 – 10 min	30 Hz
Fresh, frozen, dried, and ethanol preserved insects, e.g., <i>Drosophila melanogaster</i>	other device	to be optimized by user	see recommendations below

*Note:* Stability testing has been conducted on the NucleoSpin® Bead Tubes Type D with the MN Bead Tube Holder on a Vortex-Genie® 2 and with a mixer mill MM300 (Retsch®) at highest frequency (30 Hertz). NucleoSpin® Bead Tubes Type D withstand shaking for several hours in the MN Bead Tube Holder on a Vortex-Genie® 2 and for up to 30 minutes on a mixer mill MM300 (Retsch®) at highest frequency (30 Hertz).

For optimal sample processing, avoidance of DNA fragmentation, and highest DNA yield see table above for recommendations of adequate disruption times. Other disruption devices (see section 2.4.1) will require different settings regarding frequency and duration for optimal performance with the selected sample material. Please note that the position of the tube within the machine (mixer mill, Retsch®) is important for optimal performance! Please refer to the instrument manual of the disruption device.

**WARNING:** Many disruption devices can cause high mechanical stress on the bead tubes. Depending on bead tube type and content (beads, liquid volume, sample type), especially high frequency of shaking and / or long shaking duration can cause breaking up of the bead tubes. **It is the responsibility of the user to perform initial stability tests to ensure stability of bead tubes during the individual experimental setup!**

These tests should be performed with water instead of lysis buffer in order to avoid spillage of chaotropic lysis buffer in case of tube breakage. Integrity and tightness of the tube need to be controlled after every run.

**WARNING:** *In section 5 a certain liquid volume during disruption is recommended. The reduction of liquid content will severely increase the mechanical impact by the steel beads and can result in damage of DNA and tube.*

## 2.5 Elution procedures

In addition to the standard method, several modifications are possible to increase yield, concentration, and convenience.

- **Convenient elution (standard elution):** Elution can be performed by a single addition of 100  $\mu\text{L}$  Elution Buffer onto the column.
- **High yield:** Elution can be performed in two serial elutions of 100  $\mu\text{L}$  each, resulting in a total volume of 200  $\mu\text{L}$ .
- **High concentration:** Elution can be performed by application of 100  $\mu\text{L}$  Elution Buffer, which is then re-used in a second elution step, resulting in 100  $\mu\text{L}$  eluate with high DNA concentration. Alternatively, the elution volume can be reduced down to 25  $\mu\text{L}$ .

### 3 Storage conditions and preparation of working solutions

**Attention:**

*Lysis Buffer MG and Wash Buffer BW contain chaotropic salts! Wear gloves and goggles!*

CAUTION: Buffers MG and BW contain chaotropic salts which can form highly reactive compounds when combined with bleach (sodium hypochlorite). DO NOT add bleach or acidic solutions directly to the sample-preparation waste!

All kit components can be stored at room temperature (15–25 °C) and are stable until: see package label.

Before starting any **NucleoSpin® DNA Insect** protocol, prepare the following:

- **Wash Buffer B5:** Add the indicated volume of ethanol (96–100 %) to **Wash Buffer B5 Concentrate**. Mark the label of the bottle to indicate that ethanol was added. Wash Buffer B5 can be stored at 15–25 °C.
- **Liquid Proteinase K** is ready to use. After first time use, store Liquid Proteinase K at 4 °C or -20 °C.

NucleoSpin® DNA Insect		
REF	10 preps 740470.10	50 preps 740470.50
Wash Buffer B5 (Concentrate)	6 mL Add 24 mL ethanol	6 mL Add 24 mL ethanol

## 4 Safety instructions

When working with the **NucleoSpin® DNA Insect** kit wear suitable protective clothing (e.g., lab coat, disposable gloves, and protective goggles). For more information consult the appropriate Material Safety Data Sheets (MSDS available online at [www.mn-net.com/msds](http://www.mn-net.com/msds)).



Caution: Guanidine hydrochloride in buffer BW and guanidinium thiocyanate in buffer MG can form highly reactive compounds when combined with bleach! Thus, do not add bleach or acidic solutions directly to the sample preparation waste.

The waste generated with the **NucleoSpin® DNA Insect** kit has not been tested for residual infectious material. A contamination of the liquid waste with residual infectious material is highly unlikely due to strong denaturing lysis buffer and Proteinase K treatment but it cannot be excluded completely. Therefore, liquid waste must be considered infectious and should be handled and discarded according local safety regulations.

### 4.1 Disposal

Dispose hazardous, infectious or biologically contaminated materials in a safe and acceptable manner and in accordance with all local and regulatory requirements.

## 5 Protocols

### 5.1 Protocol for fresh, frozen, dried, and ethanol preserved insect samples

#### Before starting the preparation:

- Check if Buffer B5 was prepared according to section 3.
- Check section 2.4 for lysis and disruption of sample material.

#### 1 Prepare sample

Place the insect sample to a NucleoSpin® Bead Tube Type D (provided).



+ 100 µL BE

*Up to approx. 40 mg of wet weight insect sample can be processed. Remove excess liquid (e.g., water, ethanol) from the sample, with a filter paper.*

Add **100 µL Elution Buffer BE** to the sample.

*Alternatively, high quality grade water (not provided) can be used.*

#### 2 Lyse sample

Add **40 µL Buffer MG**.

Then, add **10 µL Liquid Proteinase K** and close the tube.



+ 40 µL MG

+ 10 µL Liquid  
Proteinase K

**Agitate** the NucleoSpin® Bead Tube in the MN Bead Tube Holder on a Vortex-Genie® 2. Alternatively a swing mill (Retsch®) can be used (see section 2.4.3).

**Agitate**

*Note: Optimal agitation duration, speed/frequency depends on the device used. For the MN Bead Tube Holder it is approximately 20 min; in a mixer mill MM200, MM300, MM400 (Retsch®), e.g., 0.5–5 min at maximal frequency (30 Hertz) is suitable (see section 2.4). On a swing mill, position of the tube in the mill can considerably influence the result. Please refer to the instruction manual of the device used. Respect warnings in section 1.2 and 2.4.3 if other devices are intended to be used!*

Centrifuge the NucleoSpin® Bead Tube **30 s** at **11,000 x g** to clean the lid.



**11,000 x g,  
30 s**

*Note:* In this step foam is removed from the screw cap to allow clean opening of the tube.

**Attention:** Do neither centrifuge at higher g-force, nor longer because this might damage the NucleoSpin® Bead Tubes.

### 3 Adjust DNA binding conditions

Add **600 µL Buffer MG** and **mix** (e.g, vortex for 3 s).



**+ 600 µL MG  
Mix**

*Note:* Steel beads should be agitated in the tube; some residual pellet (cell debris) may remain on the bottom of the tube.

Centrifuge for **30 s** at **11,000 x g**.



**11,000 x g,  
30 s**

*Note:* This centrifugation step is performed in order to clean the lid and sediment steel beads and cell debris.

**Attention:** Do neither centrifuge at higher g-force, nor longer because this might damage the NucleoSpin® Bead Tubes

### 4 Bind DNA

Transfer the supernatant (~500–600 µL) onto the NucleoSpin® DNA Insect Column, placed in a 2 mL Collection Tube (provided).



**Load samples**

Centrifuge for **30 s** at **11,000 x g**. Discard Collection Tube with flow through. Put column into a fresh Collection Tube (2 mL, provided).



**11,000 x g  
30 s**

### 5 Wash silica membrane

#### **1<sup>st</sup> wash**

Add **500 µL Buffer BW**. Centrifuge for **30 s** at **11,000 x g**. Discard flowthrough and place the column back into the Collection Tube.

**+ 500 µL BW**

**11,000 x g,  
30 s**

**+ 500 µL B5**

#### **2<sup>nd</sup> wash**

Add **500 µL Buffer B5** to the column and centrifuge for **30 s** at **11,000 x g**. Discard flowthrough and place the column back into the Collection Tube.

**11,000 x g,  
30 s**

**6 Dry silica membrane**

Centrifuge the column for **30 s** at **11,000 x g**.

*Note:* Residual wash buffer is removed in this step.



**11,000 x g,  
30 s**

**7 Elute highly pure DNA**

Place the NucleoSpin® DNA Insect Column into a 1.5 mL nuclease-free tube (not provided) and add **100 µL Elution Buffer BE** onto the column. Incubate at **room temperature** for **1 min**. Centrifuge **30 s** at **11,000 x g**.

*For alternative elution procedures see section 2.5*



**100 µL BE  
RT, 1 min**



**11,000 x g,  
30 s**

## 5.2 Protocol for purification of DNA from hard to lyse bacteria in insect samples

The purification of DNA from hard to lyse bacteria (e.g., gram positive bacteria) in insect samples can be challenging as disruption of the two organisms require individual mechanical forces.

MACHEREY-NAGEL has therefore developed the NucleoSpin® Bead Tubes Type E, which contain 40–400 µm glass beads as well as 3 mm steel beads. NucleoSpin® Bead Tubes Type E can be used according to the protocol described for NucleoSpin® Bead Tubes Type D in section 5.

However, the use of NucleoSpin® Bead Tubes Type E is a very harsh method in terms of sample disruption. Please note that processing time on a selected disruption device (e.g., MN Bead Tube Holder or mixer mill (Retsch®)) has to be optimized by the user with regard to sample type, amount, and downstream application. Long disruption duration on high impact machines (e.g., mixer mill (Retsch®)) can cause total DNA loss due to massive DNA fragmentation. Please contact the technical service for further information.

Respect warnings in section 1.2 and 2.4.3 when using NucleoSpin® Bead Tubes Type E.

## 6 Appendix

### 6.1 Troubleshooting

Problem	Possible cause and suggestions
Damaged bead tubes	<p data-bbox="330 311 988 335"><i>Unsuitable disruption device or intensity</i></p> <ul data-bbox="330 343 988 438" style="list-style-type: none"> <li data-bbox="330 343 988 438">• High force disruption devices can damage NucleoSpin® Bead Tubes Type D and E. Respect warnings in section 1.2 and 2.4.3. Use the recommended MN Bead Tube Holder.</li> </ul>
	<p data-bbox="330 454 988 478"><i>Incomplete lysis</i></p> <ul data-bbox="330 494 988 550" style="list-style-type: none"> <li data-bbox="330 494 988 550">• Adjust lysis conditions (bead tube type, agitation device, duration, or frequency).</li> </ul> <p data-bbox="330 566 988 590"><i>Reagents not applied properly</i></p> <ul data-bbox="330 606 988 630" style="list-style-type: none"> <li data-bbox="330 606 988 630">• Prepare Buffer B5 according to instructions (section 3).</li> </ul>
No or poor DNA yield	<p data-bbox="330 646 988 670"><i>Suboptimal elution of DNA from the column</i></p> <ul data-bbox="330 694 988 997" style="list-style-type: none"> <li data-bbox="330 694 988 774">• For certain sample types, preheat Buffer BE to 70 °C before elution. Apply Buffer BE directly onto the center of the silica membrane.</li> <li data-bbox="330 790 988 869">• Elution efficiencies decrease dramatically if elution is done with buffers with a pH &lt; 7.0. Use slightly alkaline elution buffers like Buffer BE (pH 8.5).</li> <li data-bbox="330 885 988 997">• Especially when expecting high yields from large amounts of material, we recommend elution with 200 µL Buffer BE and incubation of the closed columns in an incubator at 70 °C for 5 min before centrifugation.</li> </ul>
	<p data-bbox="330 1013 988 1037"><i>High amount of sample material</i></p> <ul data-bbox="330 1053 988 1276" style="list-style-type: none"> <li data-bbox="330 1053 988 1276">• If using more than 20 mg of sample material, DNA yield can be increased by addition of 20 µL ethanol after the addition of 600 µL Buffer MG in step 3. However, ethanol addition slightly increases RNA copurification. In order to avoid RNA copurification due to ethanol addition, incubate lysate after disruption for 5 min at 70 °C, then add 600 µL Buffer MG and 20 µL ethanol, mix, and proceed with the transfer onto the NucleoSpin® DNA Insect Column.</li> </ul>

Problem	Possible cause and suggestions
Poor DNA quality	<p data-bbox="333 209 519 228"><i>High <math>A_{260}/A_{280}</math> ratio</i></p> <ul data-bbox="333 252 978 520" style="list-style-type: none"> <li data-bbox="333 252 978 520">• Ratios &gt; 1.9 can be caused by RNA contamination. Usually, such RNA contamination does not interfere with downstream application. Depending on sample type, amount, and disruption procedure, preparations might contain small amounts of RNA. If it is necessary to reduce RNA contamination to the lowest possible level, incubate the lysate after the disruption step for 5 min at 70 °C in order to inactivate the Proteinase K. After cooling to room temperature, add 20 µL RNase A (20 mg/mL) and incubate 5 min. Continue with the application of the lysate onto the column.</li> </ul> <p data-bbox="333 544 628 563"><i>Reagents not applied properly</i></p> <ul data-bbox="333 587 945 606" style="list-style-type: none"> <li data-bbox="333 587 945 606">• Prepare Buffer B5 according to instructions (see section 3).</li> </ul>
Clogged columns	<p data-bbox="333 632 645 651"><i>Too much sample material used</i></p> <ul data-bbox="333 675 978 748" style="list-style-type: none"> <li data-bbox="333 675 978 748">• Make sure to centrifuge the lysate after cell disruption in order to sediment beads and cell debris. Only transfer cleared supernatant onto the column.</li> </ul>
Suboptimal performance of genomic DNA in enzymatic reactions	<p data-bbox="333 775 605 794"><i>Carry-over of ethanol or salt</i></p> <ul data-bbox="333 818 978 1018" style="list-style-type: none"> <li data-bbox="333 818 978 922">• Make sure to centrifuge <math>\geq 1</math> min at 11,000 x g in order to remove all of ethanolic Buffer B5 before eluting the DNA. If, for any reason, the level of Buffer B5 has reached the column outlet after drying, repeat the centrifugation.</li> <li data-bbox="333 946 978 1018">• Do not chill Buffer B5 before use. Cold buffer will not remove salt effectively. Equilibrate Buffer B5 to room temperature before use.</li> </ul> <p data-bbox="333 1042 810 1061"><i>Contamination of DNA with inhibitory substances</i></p> <ul data-bbox="333 1085 978 1131" style="list-style-type: none"> <li data-bbox="333 1085 978 1131">• Do not elute DNA with TE buffer. EDTA may inhibit enzymatic reactions. Repurify DNA and elute in Buffer BE.</li> </ul>

## 6.2 Ordering information

Product	REF	Pack of
NucleoSpin® DNA Insect	740470.10 / .50	10/50 preps
MN Bead Tube Holder	740469	1 piece
NucleoSpin® Soil	740780.10 / .50 / .250	10/50/250 preps
NucleoSpin® DNA Lipid Tissue	740471.10 / .50	10/50 preps
NucleoSpin® Microbial DNA	740235.10 / .50	10/50 preps
NucleoSpin® Bead Tubes Type A (0.6–0.8 mm ceramic beads, recommended for soil and sediments)	740786.50	50 pieces
NucleoSpin® Bead Tubes Type B (40–400 µm glass beads, recommended for bacteria)	740812.50	50 pieces
NucleoSpin® Bead Tubes Type C (1–3 mm corundum, recommended for yeast)	740813.50	50 pieces
NucleoSpin® Bead Tubes Type D (3 mm steel beads, recommended for insects)	740814.50	50 pieces
NucleoSpin® Bead Tubes Type E (40–400 µm glass beads and 3 mm steel beads, recommended for hard to lyse bacteria within insect samples)	740815.50	50 pieces
NucleoSpin® Bead Tubes Type F (1–3 mm corundum and 3 mm steel beads, use only with MN Bead Tube Holder!)	740816.50	50 pieces
NucleoSpin® Bead Tubes Type G (5 mm steel beads, recommended for plant material)	740817.50	50 pieces
Buffer BE	740306.100	125 mL
Buffer B5 Concentrate (for 125 mL Buffer B5)	740921	25 mL
Buffer BW	740922	100 mL

<b>Product</b>	<b>REF</b>	<b>Pack of</b>
Liquid Proteinase K	740396	5 mL
RNase A	740505.50 740505	50 mg 100 mg
Collection Tubes (2 mL)	740600	1000

Visit [www.mn-net.com](http://www.mn-net.com) for more detailed product information.

### **6.3 Product use restriction/warranty**

All MACHEREY-NAGEL products are designed for their intended use only. They are not intended to be used for any other purpose. The description of the intended use of the products can be found in the original MACHEREY-NAGEL product leaflets. Before using our products, please observe the instructions for use and the safety instructions from the respective Material Safety Data Sheet of the product.

This MACHEREY-NAGEL product is carrying documentation stating specifications and other technical information. MACHEREY-NAGEL warrants to meet the stated specifications. The provided warranty is limited to the data specifications and descriptions as given in the original MACHEREY-NAGEL literature. No other statements or representations, written or oral, by MACHEREY-NAGEL's employees, agents or representatives, except written statements signed by a duly authorized officer of MACHEREY-NAGEL are authorized. They should not be relied upon by the customer and are not a part of a contract of sale or of this warranty.

Liability for all possible damages that occur in any connection with our products is limited to the utmost minimum as stated in the general business terms and conditions of MACHEREY-NAGEL in their latest edition which can be taken from the company's website. MACHEREY-NAGEL does not assume any further warranty.

Products and their application are subject to change. Therefore, please contact our Technical Service Team for the latest information on MACHEREY-NAGEL products. You may also contact your local distributor for general scientific information. Descriptions in MACHEREY-NAGEL literature are provided for informational purposes only.

Last updated: 08/2022, Rev. 04

Please contact:  
MACHEREY-NAGEL GmbH & Co. KG  
Tel.: +49 24 21 969-333  
[support@mn-net.com](mailto:support@mn-net.com)


---

Trademarks:

Bullet Blender is a registered trademark of Next Advance  
FastPrep is a registered trademark of MP Biomedicals, LLC  
MagNA Lyser is a trademark of Roche Diagnostics GmbH  
Mini-Beadbeater is a trademark of Biospec Products  
NucleoSpin is a registered trademark of MACHEREY-NAGEL GmbH & Co KG  
Precellys is a registered trademark of Bertin Technologies  
Retsch is a registered trademark of Retsch GmbH  
Vortex-Genie is a registered trademark of Scientific Industries

All used names and denotations can be brands, trademarks, or registered labels of their respective owner – also if they are not special denotation. To mention products and brands is only a kind of information (i.e., it does not offend against trademarks and brands and can not be seen as a kind of recommendation or assessment). Regarding these products or services we can not grant any guarantees regarding selection, efficiency, or operation.



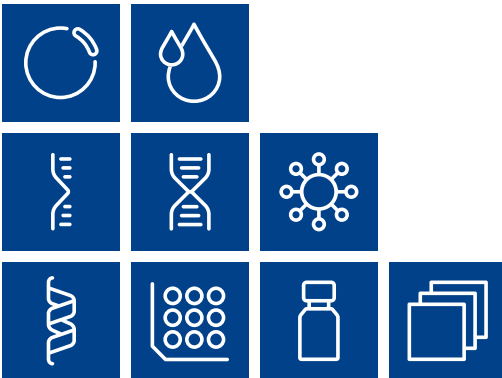

Plasmid DNA

Clean up

RNA

DNA


Viral RNA and DNA

Protein

High throughput

Accessories

Auxiliary tools

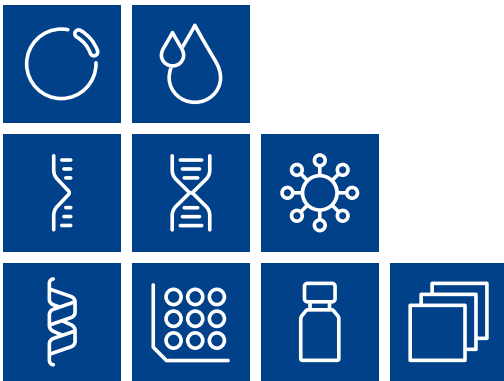

# MACHEREY-NAGEL

[www.mn-net.com](http://www.mn-net.com)

MACHEREY-NAGEL GmbH & Co. KG · Valencienner Str. 11 · 52355 Düren · Germany

DE +49 24 21 969-0 [info@mn-net.com](mailto:info@mn-net.com)

CH +41 62 388 55 00 [sales-ch@mn-net.com](mailto:sales-ch@mn-net.com)

FR +33 388 68 22 68 [sales-fr@mn-net.com](mailto:sales-fr@mn-net.com)

US +1 888 321 62 24 [sales-us@mn-net.com](mailto:sales-us@mn-net.com)

