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Description:

pH-Fix 0–14 are test strips to determine the pH of aqueous samples in the range of pH 0–14.

They are suitable for visual evaluation using a color scale and for reflectometric evaluation with the QUANTOFIX® Relax.

The indicator is covalently bound to the test fields and does not bleed into the sample.

Pack content:

- 100 test strips (6 × 85 mm)
- 1 reference color scale

Measuring range:

Visually:

pH - 0 · 1 · 2 · 3 · 4 · 5 · 6 · 7 · 8 · 9 · 10 · 11 · 12 · 13 · 14

QUANTOFIX® Relax:

pH - 1–13

Additional indications:

The base color of the paper in its dry state has no effect on functionality.

Hazard warnings:

This test does not contain hazardous substances that must be labelled.

Instructions for use (visual):

- 1 Dip the test strip with all test fields into the sample for 1 sec.
- 2 Shake off excess liquid.
- 3 Compare the test field to the color scale.

In cases of insufficient buffering, leave the test strip in the sample until the colors remain constant.

Instructions for use (QUANTOFIX® Relax):

- 1 Select the [pH 0–14] method.
- 2 Dip the test strip with all test fields into the sample for 1 sec.
- 3 Shake off excess liquid.
- 4 Place the test strip on the test strip holder.
- 5 Read the test results once the reaction time has elapsed.

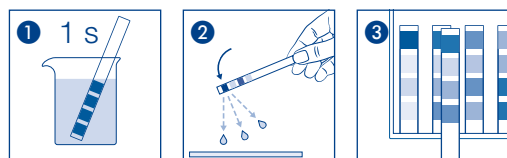
Disposal:

Dispose of used test strips in household waste.

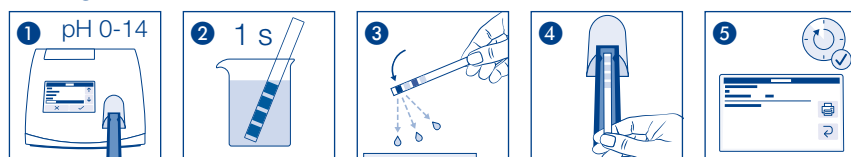
Storage:

Protect test strips from sunlight and moisture. Store the container in a cool, dry place (storage temperature 4–30 °C / 39–86 °F). If stored properly, the test strips are stable until the printed expiration date.

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MACHEREY-NAGEL GmbH & Co. KG · Valencienner Str. 11 · 52355 Düren · Germany