

Customer information

Method update NANOCOLOR® UV/VIS II and VIS II

Overview

Photometer	NANOCOLOR® UV/VIS II und VIS II
Method number	r1382
Availability	from 10 th of September 2026
Price	free of charge
Required accessories	USB stick




Important note

Even if an update does not influence the data and special methods in the device, MACHEREY-NAGEL recommends some security actions to prevent data loss. Before performing an update save the memory and the database of the device via the -icon in the corresponding menu. Furthermore, a regular backup via the backup menu in the settings is recommended. It is recommended to save special methods on a USB stick prior to the update process. A more detailed description of the data saving is provided in the manual of the device.



Modifications

- Optimization of Method 0-10 Vicinal Diketones on the NANOCOLOR® UV/VIS II



Update procedure

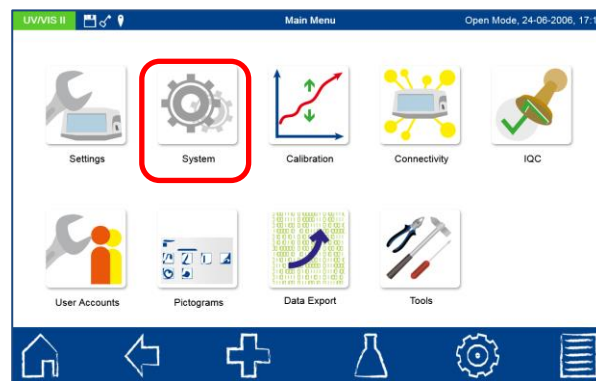

- Copy the method update file to the root directory of a USB flash drive and connect the USB drive to one of the photometer's USB-A ports. Then navigate to the "System" menu via the main menu

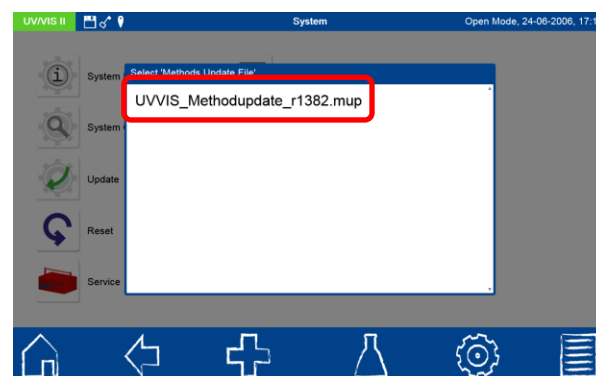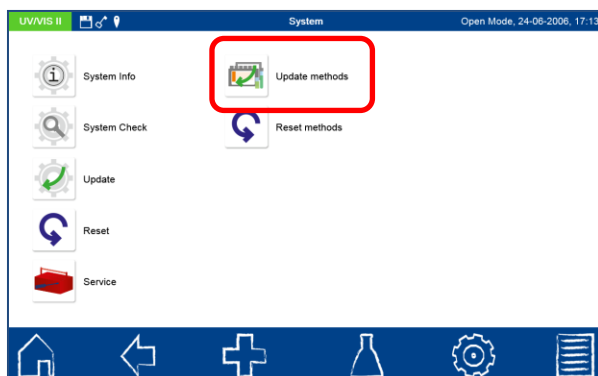

- Open the "Method Update" menu and select the desired method update file (.mup file)

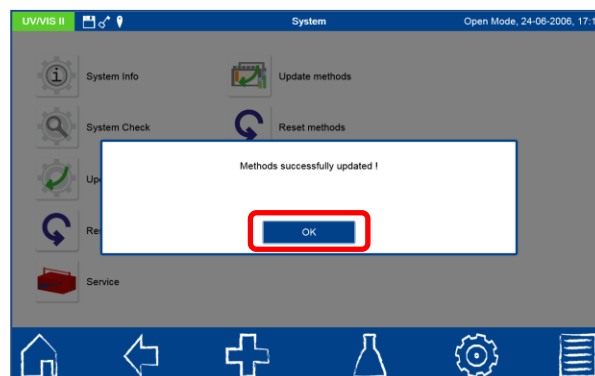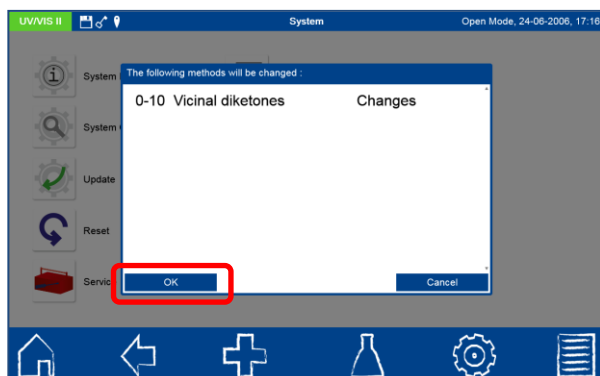

- On the following screen, all methods that will be modified or added by the update are displayed. Confirm with “OK” to perform the update, or select “Cancel” to abort the process.
- Once the update has been completed successfully, which may take several minutes, confirm the message

Contact

If you have further questions, please do not hesitate to contact us:

Technical support and customer service

Tel.: +49 24 21 969-333
E-Mail: support@mn-net.com

For orders (international)

E-Mail: sales@mn-net.com