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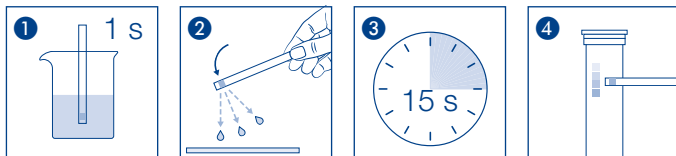
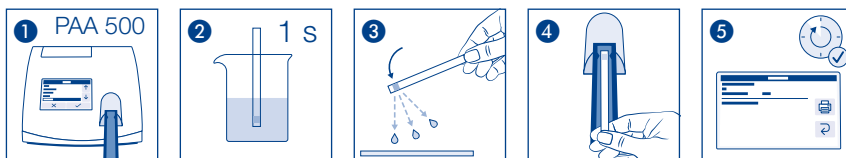
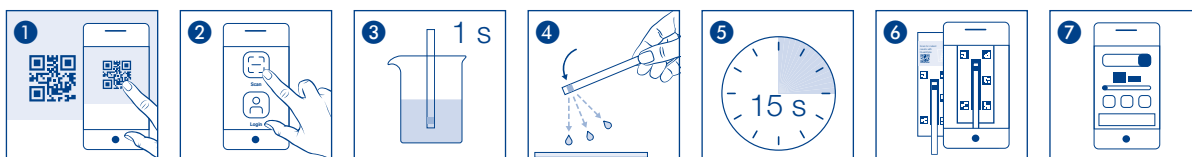
Rev. 10/2026

Description:

QUANTOFIX® Peracetic acid 500 are test strips for the semi-quantitative determination of peracetic acid (PAA) (CH_3COOOH) in solutions.

They are suitable for visual evaluation using a color scale, for reflectometric evaluation with the QUANTOFIX® Relax, and for evaluation using the *QuantifyGo*® web application on a smartphone.

Pack content:

- 1 aluminium container with 100 test strips and reference color scale

Measuring range:

Visually:

Peracetic acid - 0 · 50 · 100 · 200 · 300 · 400 · 500 mg/L PAA

QUANTOFIX® Relax:

Peracetic acid - 50 – 500 mg/L PAA

QuantifyGo®:

Peracetic acid - 50 – 500 mg/L PAA

Method:

Peracetic acid - Redox reaction

Color change:

Peracetic acid - yellow → green

General indications:

Always remove only the number of test strips needed. Immediately close the container after removing test strips. Do not touch the test field.

Hazard warnings:

This test does not contain hazardous substances that must be labelled.

Instructions for use (visual):

- 1 Dip the test strip into the sample for 1 sec.
- 2 Shake off excess liquid.
- 3 Wait 15 sec.
- 4 Then compare the test field immediately with the color scale and read the value that most closely matches the color of the test field (reading accuracy: $\pm \frac{1}{2}$ color field)

For interpreting the results, do not take color changes into consideration after the reaction time has elapsed.

Instructions for use (QUANTOFIX® Relax):

- 1 Select the [PAA 500] method.
- 2 Dip the test strip into the sample for 1 sec.
- 3 Shake off excess liquid.
- 4 Place the test strip on the test strip holder.
- 5 Read the test results once the reaction time has elapsed.

Instructions for use (*QuantifyGo*®):

- 1 Visit <https://quantifygo.mn-net.com/FNZXKX> or scan the QR code.
- 2 Press the "Scan" button in the *QuantifyGo*® application.
- 3 Dip the test strip into the sample for 1 sec.
- 4 Shake off excess liquid.
- 5 Wait 15 sec.
- 6 Place the test strip on the marked area and take a photo.
- 7 The result appears immediately on the screen.

Quality control:

To check the correct functioning of the test strips, use a peracetic acid solution with a concentration of 200 mg/L. To prepare the solution, use a volumetric flask and add to 0.11 mL peracetic acid 39 % fully distilled water, filling the flask to the 250 mL mark. Immediately perform the measurement with the test strip. If the control solution produces a negative result even after repeating the process, then the remaining unused test strips must be discarded. Even during a negative control (inserting a test strip into fully distilled water), no color deviation from the 0 mg/L PAA test field may occur. Possible reasons for incorrect functioning of the test strips may be that the use-by-date has been exceeded, the container has been left open for too long or has been stored incorrectly.

Interferences:

If the sample solution has a pH value of 4–6, the reaction will take place without interferences. Acidic solutions must be buffered with sodium acetate, and alkaline solutions with diluted acetic acid to a pH of 5. The following ions interfere with the determination when the concentrations below are exceeded (tested with 0 and 200 mg/L respectively of peracetic acid solutions): 5 mg/L: bound chlorine (chloramine), Fe^{2+} 10 mg/L: free chlorine (hypochlorite) 15 mg/L: ascorbic acid, Fe^{3+} 75 mg/L: sulphite (SO_3^{2-}) 1000 mg/L: formaldehyde, nitrate (NO_3^-), hydrogen peroxide (H_2O_2) 30 °d: total hardness

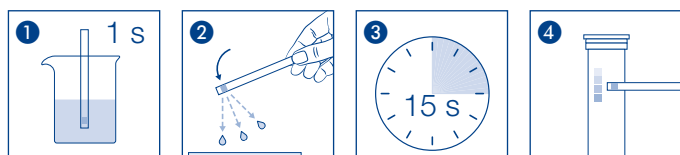
Disposal:

Dispose of used test strips in household waste.

Storage:

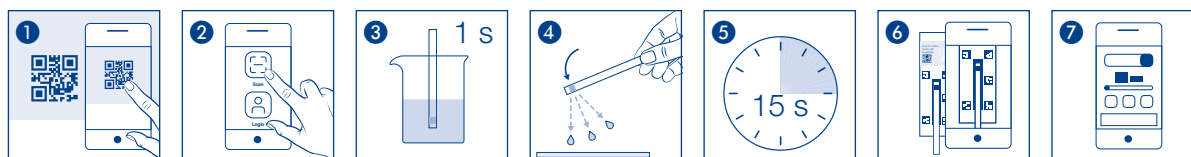
Protect test strips from sunlight and moisture. Store the container in a cool, dry place (storage temperature 4–30 °C / 39–86 °F). If stored properly, the test strips are stable until the printed expiration date.

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MACHEREY-NAGEL GmbH & Co. KG · Valencienner Str. 11 · 52355 Düren · Germany