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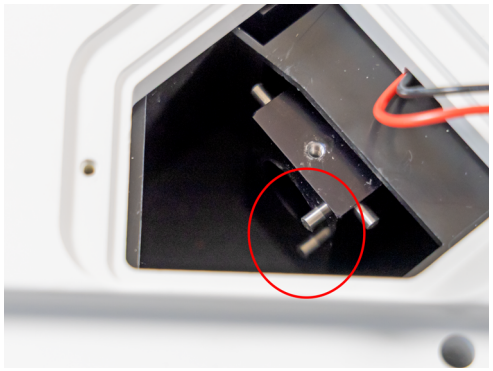
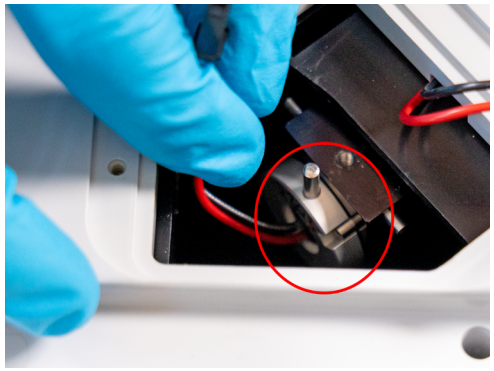
1. Important information

A halogen lamp provides the light for measurements in the 320 nm – 1100 nm range. For each standard measurement, the lamp emits a short light pulse. This results in very low energy consumption and a long lamp life cycle.

DANGER

Risk of electric shock: Always disconnect the device from the power source before replacing the lamp. Do not reconnect the device to the power source until the lamp replacement is complete.

WARNING

Risk of burns. The light source becomes very hot when the device is in use. There is a risk of burns. You should therefore allow the lamps to cool down for at least 30 minutes before replacing them.

2. Required materials

- Screwdriver
- Soft, safe work surface
- Halogen lamp REF 919759

NOTE


Use only accessories and replacement parts approved by the manufacturer. If parts not approved by the manufacturer are used, the operator is responsible for the device's compliance, and the warranty will be void

3. Implementation

3.1 Unplug the device from the power source and place it upside down on a soft work surface.

3.2 Remove the screw from the lamp cover and open it.

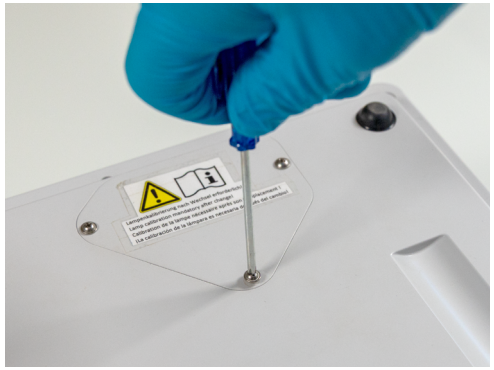MACHEREY-NAGEL

www.mn-net.com


MACHEREY-NAGEL GmbH & Co. KG · Valencienner Str. 11 · 52355 Düren · Germany

DE +49 24 21 969-0 info@mn-net.com
CH +41 62 388 55 00 sales-ch@mn-net.com
FR +33 388 68 22 68 sales-fr@mn-net.com
US +1 888 321 62 24 sales-us@mn-net.com

Replacing the halogen lamp for the *NANOCOLOR*® Advance spectrophotometer

3.3 Remove the screw from the metal cover above the halogen bulb and take the cover off.

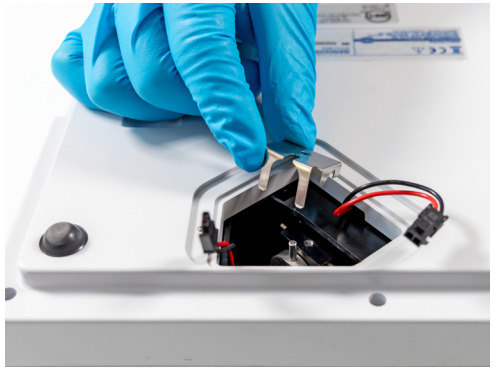
3.4 Unplug the halogen lamp from the outlet and remove it from its holder.

3.5 Insert a new halogen bulb using the guide pins and connect the plug.

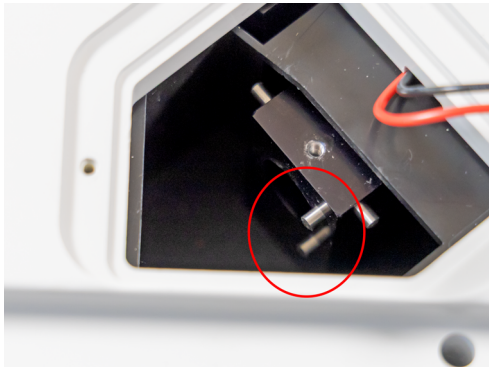3.6 Secure the metal trim and the cover using the screws.

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