

NucleoMag[®] RNA Pro– NucleoMag[®] X32

Protocol details

Application	Highly pure RNA from tissue and plant leaf material
Kit	NucleoMag [®] RNA Pro
REF	744360 (.1 / .4)
Protocol name	NMRNAProR1
Run time	approx. 77 min



Protocol steps

Procedure	
1	Perform lysis according to the user manual NucleoMag [®] RNA Pro.
2	Fill the 96-well Deep-well plates according to the table below.
3	Load the plates on the NucleoMag [®] X32.
4	Insert tip combs on the mounting grooves.*
5	Select the protocol from the instrument menu and start protocol NMRNAProR1.
6	Remove plate from the instrument when prompted on the NucleoMag [®] X32.
7	Dispense 350 µL of Binding Buffer MRB to the DNA digestion (column 4 + 10).
8	Place the plate back into the NucleoMag [®] X32 and proceed with the protocol

Notes: * Please equip all tip combs in order to cover the magnetic rods in used and unused wells.
The protocol includes a rDNase incubation step that requires a manual intervention (addition of Binding Buffer MRB)

Required consumables and instrumentation provided by user

Product	Specification	REF
NucleoMag X32 [®]	Automated nucleic extraction system for MACHEREY-NAGEL's NucleoMag [®] kits enabling parallel processing of up to 32 samples	747020
96 Deep-well Plate	96 deep-well plates for NucleoMag [®] X32 (25 pieces)	744955
Tip Combs	8-well tip combs for NucleoMag [®] X32 (50 pieces)	744960

Find more NucleoMag[®] X32 protocol information and scripts here: <https://www.mn-net.com/de/NucleoMag-X32-Scripts>



Loading table


Position	Reagents	Samples per plate
Column 1 + 7	Cleared Lysate (350 µL), Binding Reagent (250 µL)*, NucleoMag® B-Beads (20 µL)	Sample 1-8 Sample 9-16
Column 2 + 8	Wash Buffer MRW (900 µL)	Sample 1-8 Sample 9-16
Column 3 + 9	Ethanol 70% (900 µL)	Sample 1-8 Sample 9-16
Column 4 + 10	rDNase reaction mixture (300 µL)	Sample 1-8 Sample 9-16
Column 5 + 11	Ethanol 70% (900 µL)	Sample 1-8 Sample 9-16
Column 6 + 12	Elution Buffer MRE (100 µL)	Sample 1-8 Sample 9-16

Notes:

Please refer to the image below for a visual representation of the loading scheme

* Use binding reagent for different sample types according to the user manual NucleoMag® RNA Pro.

Loading scheme

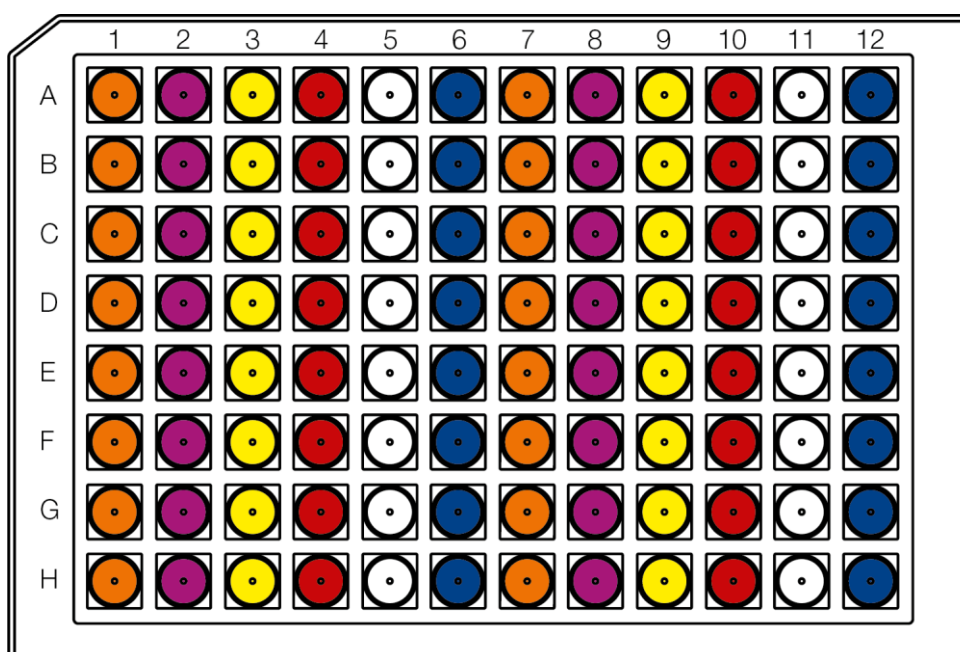

Disclaimer

Information

MACHEREY-NAGEL GmbH & Co. KG makes every effort to include accurate and up-to-date information within this publication; however, it is possible that omissions or errors might have occurred. MACHEREY-NAGEL GmbH & Co. KG cannot, therefore, make any representations or warranties, expressed or implied, as to the accuracy or completeness of the information provided in this publication. Changes in this publication can be made at any time without notice. For technical details and detailed procedures of the specifications provided in this document please contact your MACHEREY-NAGEL representative. This publication may contain reference to applications and products which are not available in all markets. Please check with your local sales representative. All mentioned trademarks are protected by law. All used names and denotations can be brands, trademarks, or registered labels of their respective owner – also if they are not special denotation. To mention products and brands is only a kind of information (i.e., it does not offend against trademarks and brands and can not be seen as a kind of recommendation or assessment). Reference herein to any specific commercial product, process, or service by trade name, trademark, manufacturer, or otherwise, does not necessarily constitute or imply its endorsement, recommendation, or support by MACHEREY-NAGEL GmbH & Co. KG. Any views or opinions expressed herein by the authors' do not necessarily state or reflect those of MACHEREY-NAGEL. NucleoMag® is a registered trademark of MACHEREY-NAGEL GmbH & Co. KG, Düren, Germany.