



visocolor® ECO Aluminium

Aluminum / Aluminium / Aluminio

Messbereich / Range / Domaine de mesure / Rango
0.10 – 0.50 mg/L Al³⁺

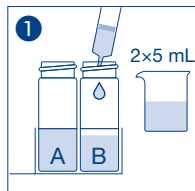


Probe
Sample
Echantillon
Muestra

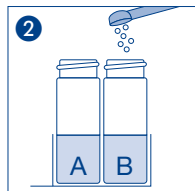


1 — Al-1




Schütteln
Shake
Agiter
Agitar

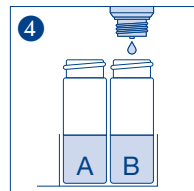

4 — Al-2

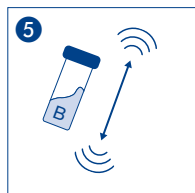


Schütteln
Shake
Agiter
Agitar

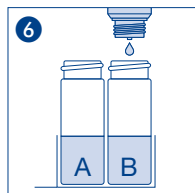


10 — Al-3



Schütteln
Shake
Agiter
Agitar

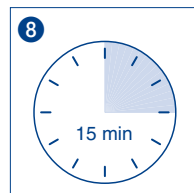

Warten
Wait
Attendre
Espere

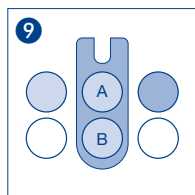

Messung
Measurement
Mesure
Medición

Geeignet für Meerwasser nach Verdünnung 1 + 9
Suitable for sea water after dilution 1 + 9
Convient à l'eau de mer après dilution 1 + 9
Adecuado para el agua de mar tras su dilución 1 + 9