

NucleoMag® Blood 200 µL – NucleoMag® X32

Protocol details

Application	DNA purification from blood
Kit	NucleoMag® Blood 200 µL
REF	744501 (.1 / .4)
Protocol name	NMBlood200
Run time	approx. 21 min



Protocol steps

Procedure	
1	Perform lysis according to the user manual NucleoMag® Blood 200 µL.
2	Fill the 96-well Deep-well plates according to the table below.
3	Load the plates on the NucleoMag® X32.
4	Insert tip combs on the mounting grooves.*
5	Select the protocol from the instrument menu and start run.

Notes: * Please equip all tip combs in order to cover the magnetic rods in used and unused wells.

Required consumables and instrumentation provided by user

Product	Specification	REF
NucleoMag X32®	Automated nucleic extraction system for MACHEREY-NAGEL's NucleoMag® kits enabling parallel processing of up to 32 samples	747020
96 Deep-well Plate	96 deep-well plates for NucleoMag® X32 (25 pieces)	744955
Tip Combs	8-well tip combs for NucleoMag® X32 (50 pieces)	744960

Find more NucleoMag® X32 protocol information and scripts here: <https://www.mn-net.com/de/NucleoMag-X32-Scripts>




Loading table

Position	Reagents	Samples per plate
Column 1 + 7	Lysate (300 µL), NucleoMag® B-Beads (25 µL), Binding Buffer MBL2 (300 µL)	Sample 1-8 Sample 9-16
Column 2 + 8	Wash Buffer MBL3 (800 µL)	Sample 1-8 Sample 9-16
Column 3 + 9	Wash Buffer MBL3 (800 µL)	Sample 1-8 Sample 9-16
Column 4 + 10	Ethanol 80% (800 µL)	Sample 1-8 Sample 9-16
Column 5 + 11	Wash Buffer MBL4 (900 µL)	Sample 1-8 Sample 9-16
Column 6 + 12	Elution Buffer MBL5 (100 µL)	Sample 1-8 Sample 9-16

Notes:

Please refer to the image below for a visual representation of the loading scheme

Loading scheme

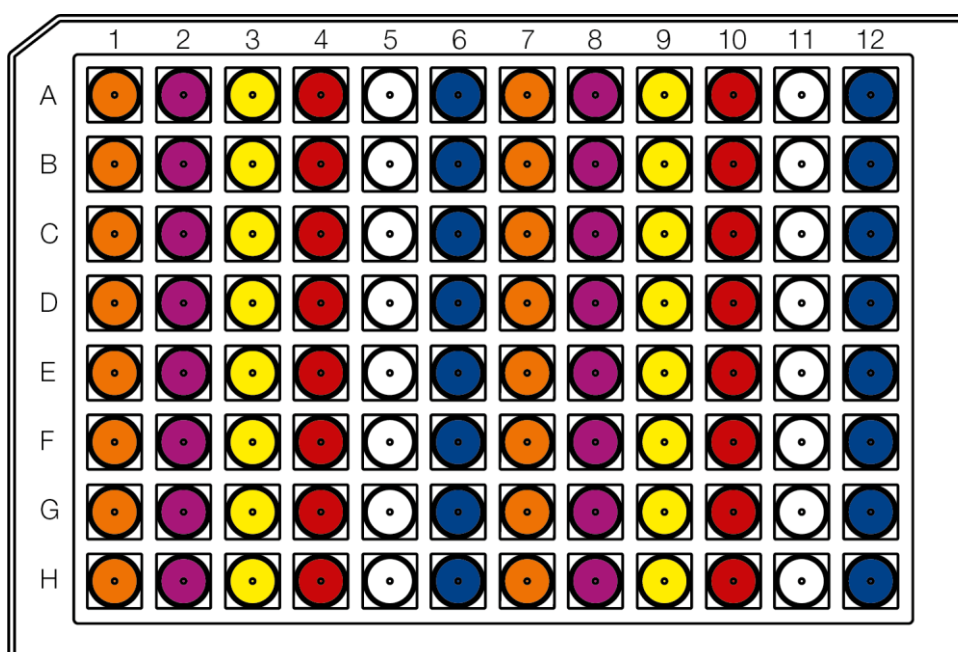

Disclaimer

Information

MACHEREY-NAGEL GmbH & Co. KG makes every effort to include accurate and up-to-date information within this publication; however, it is possible that omissions or errors might have occurred. MACHEREY-NAGEL GmbH & Co. KG cannot, therefore, make any representations or warranties, expressed or implied, as to the accuracy or completeness of the information provided in this publication. Changes in this publication can be made at any time without notice. For technical details and detailed procedures of the specifications provided in this document please contact your MACHEREY-NAGEL representative. This publication may contain reference to applications and products which are not available in all markets. Please check with your local sales representative. All mentioned trademarks are protected by law. All used names and denotations can be brands, trademarks, or registered labels of their respective owner – also if they are not special denotation. To mention products and brands is only a kind of information (i.e., it does not offend against trademarks and brands and can not be seen as a kind of recommendation or assessment). Reference herein to any specific commercial product, process, or service by trade name, trademark, manufacturer, or otherwise, does not necessarily constitute or imply its endorsement, recommendation, or support by MACHEREY-NAGEL GmbH & Co. KG. Any views or opinions expressed herein by the authors' do not necessarily state or reflect those of MACHEREY-NAGEL. NucleoMag® is a registered trademark of MACHEREY-NAGEL GmbH & Co. KG, Düren, Germany.