

NucleoMag[®] Tissue – NucleoMag[®] X32

Protocol details

Application	DNA purification from cells and tissue
Kit	NucleoMag [®] Tissue
REF	744300 (.1 / .4)
Protocol name	NMTissue
Run time	approx. 17 min



Protocol steps

Procedure	
1	Perform lysis according to the user manual NucleoMag [®] Tissue.
2	Fill the 96-well Deep-well plates according to the table below.
3	Load the plates on the NucleoMag [®] X32.
4	Insert tip combs on the mounting grooves.*
5	Select the protocol from the instrument menu and start run.

Notes: * Please equip all tip combs in order to cover the magnetic rods in used and unused wells.

Required consumables and instrumentation provided by user

Product	Specification	REF
NucleoMag X32 [®]	Automated nucleic extraction system for MACHEREY-NAGEL's NucleoMag [®] kits enabling parallel processing of up to 32 samples	747020
96 Deep-well Plate	96 deep-well plates for NucleoMag [®] X32 (25 pieces)	744955
Tip Combs	8-well tip combs for NucleoMag [®] X32 (50 pieces)	744960

Find more NucleoMag[®] X32 protocol information and scripts here: <https://www.mn-net.com/de/NucleoMag-X32-Scripts>




Loading table

Position	Reagents	Samples per plate
Column 1 + 7	Cleared lysates (225 µL), NucleoMag® B-Beads (24 µL), Binding Buffer MB2 (360 µL)	Sample 1-8 Sample 9-16
Column 2 + 8	Wash Buffer MB3 (600 µL)	Sample 1-8 Sample 9-16
Column 3 + 9	Wash Buffer MB4 (600 µL)	Sample 1-8 Sample 9-16
Column 4 + 10	Wash Buffer MB5 (900 µL)	Sample 1-8 Sample 9-16
Column 5 + 11	empty	- -
Column 6 + 12	Elution Buffer MB6 (100 µL)	Sample 1-8 Sample 9-16

Notes:

Please refer to the image below for a visual representation of the loading scheme

Loading scheme

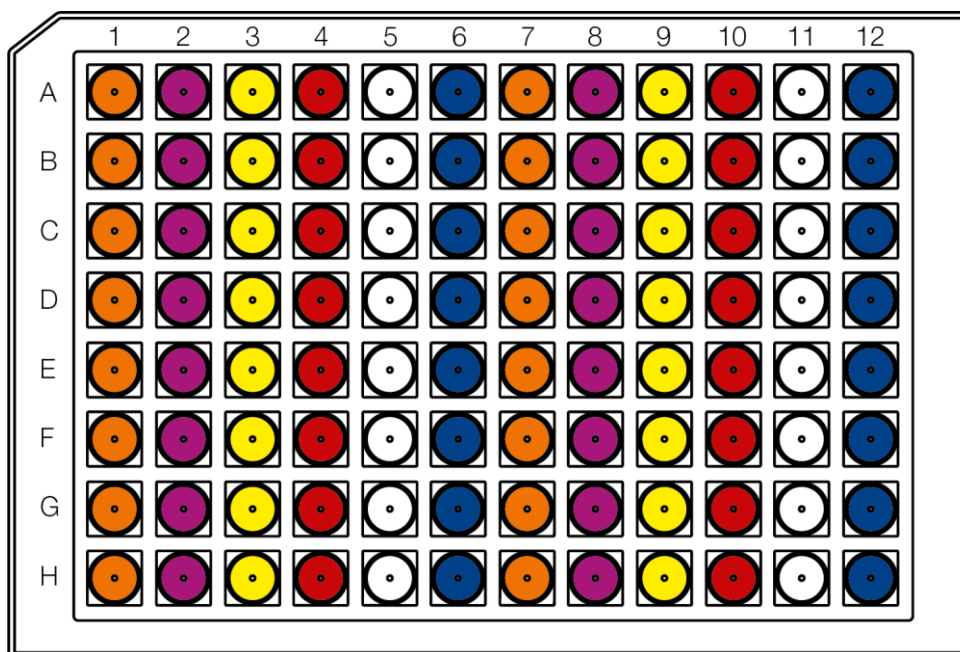

Disclaimer

Information

MACHERY-NAGEL GmbH & Co. KG makes every effort to include accurate and up-to-date information within this publication; however, it is possible that omissions or errors might have occurred. MACHERY-NAGEL GmbH & Co. KG cannot, therefore, make any representations or warranties, expressed or implied, as to the accuracy or completeness of the information provided in this publication. Changes in this publication can be made at any time without notice. For technical details and detailed procedures of the specifications provided in this document please contact your MACHERY-NAGEL representative. This publication may contain reference to applications and products which are not available in all markets. Please check with your local sales representative. All mentioned trademarks are protected by law. All used names and denotations can be brands, trademarks, or registered labels of their respective owner – also if they are not special denotation. To mention products and brands is only a kind of information (i.e., it does not offend against trademarks and brands and can not be seen as a kind of recommendation or assessment). Reference herein to any specific commercial product, process, or service by trade name, trademark, manufacturer, or otherwise, does not necessarily constitute or imply its endorsement, recommendation, or support by MACHERY-NAGEL GmbH & Co. KG. Any views or opinions expressed herein by the authors' do not necessarily state or reflect those of MACHERY-NAGEL. NucleoMag® is a registered trademark of MACHERY-NAGEL GmbH & Co. KG, Düren, Germany.