

MACHEREY-NAGEL

Sipper module

NANOCOLOR[®] FP-200

Water Analysis



Maximum performance for your lab

- Fast – Efficient work without changing the cuvette
- Precise – Increase your measuring accuracy
- Comfortable – Perfect interaction with the spectrophotometers *NANOCOLOR*[®] VIS II and *NANOCOLOR*[®] UV/VIS II

MACHEREY-NAGEL


www.mn-net.com



Sipper module NANOCOLOR® FP-200

Sipper pump for photometric analysis

Dealing with large sample numbers in photometric analysis often requires a lot of time and resources due to the preparation of the samples in different cuvettes and the frequent changes of the cell in the photometers. With the help of the sipper module NANOCOLOR® FP-200 and a corresponding flow cell (2 mm, 10 mm or 50 mm) the effort can be reduced significantly and comfortably when a high quantity of samples is present.



Using the sipper module, the danger of direct contact with chemicals is drastically reduced. A recasting of samples into the cuvettes is no longer necessary.

Perfect interaction between sipper and photometer

- Direct operation of the sipper module via the spectrophotometers NANOCOLOR® VIS II and NANOCOLOR® UV/VIS II
- Method-specific programming – no change of sipper settings necessary

Highest flexibility

- Individually adjustable pump-, rinse- and settling times
- Comfortable usage also for customized applications

Powerful analytics


- Increase of analytical speed
- Time saving for high sample throughput, especially for usage of NANOCOLOR® standard test kits

Robust system

- High chemical resistance due to special tubes
- No contact of chemicals with the pump

Ordering information

Description	REF
Sipper module NANOCOLOR® FP-200 for NANOCOLOR® UV/VIS II and NANOCOLOR® VIS II incl. manual, USB cable, tygon tube set with adapters, intake needle and certificate	919180
Tube set for sipper module NANOCOLOR® FP-200, consisting of tygon tube set with adapters and pump tubing with pre-installed fittings	919181
Manual for sipper module NANOCOLOR® FP-200	919182
Intake needle for NANOCOLOR® FP-200	919142
Stand for NANOCOLOR® FP-200	919143

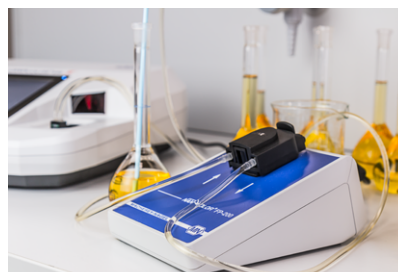

Delivery content REF 919180

www.mn-net.com

MACHEREY-NAGEL



MACHEREY-NAGEL GmbH & Co. KG
Valenciennner Str. 11
52355 Düren · Germany

DE Tel.: +49 24 21 969-0 info@mn-net.com
CH Tel.: +41 62 388 55 00 sales-ch@mn-net.com
FR Tel.: +33 388 68 22 68 sales-fr@mn-net.com
US Tel.: +1 888 321 62 24 sales-us@mn-net.com

