


visicolor[®] ECO Kieselsäure HR 200


Silica / Silice / Silice

| Method(e) / Método | | PF-3 Pool/ Drinking Water 450 nm | PF-12 ^{Plus} 450 nm | Advance 450 nm |
|--------------------|------------------------------|--|---------------------------------|-------------------|
| 5341 | 5–100 mg/L Si | ■ | ■ | ■ |
| 5342 | 10–200 mg/L SiO ₂ | ■ | ■ | ■ |

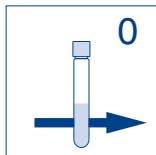


Null mit Probe
Blank with sample
Blanc avec échantillon
En blanco con muestra


Geeignet für Meerwasser
Suitable for seawater
Convient à l'eau de mer
Adecuado para el agua de mar

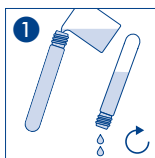


Küvetten spülen
Rinse test tube
Rincer les éprouvettes
Lavar los tubos

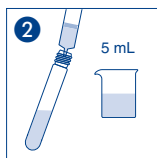


Probe
Sample
Echantillon
Muestra

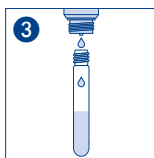


5 Δ SiO₂-1

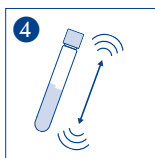

Schütteln
Shake
Agiter
Agitar

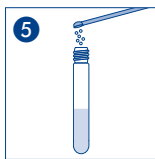


1 Δ SiO₂-2

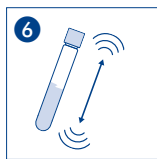


Schütteln
Shake
Agiter
Agitar




Warten
Wait
Attendre
Espere

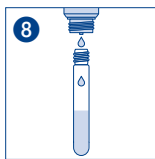

5 Δ SiO₂-3

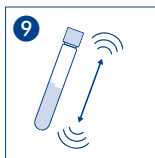


Schütteln
Shake
Agiter
Agitar

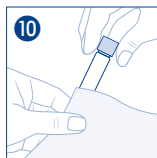

Säubern
Clean
Nettoyer
Limpiar



Warten
Wait
Attendre
Espere

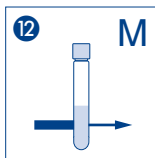

Messung
Measurement
Mesure
Medición