



visocolor® ECO Ammonium 3

Ammonium 3 / Ammonium 3 / Amonio 3

Messbereich / Range / Domaine de mesure / Rango
0.2 – 3 mg/L NH_4^+

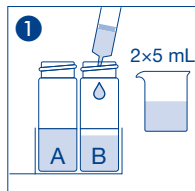


Probe
Sample
Echantillon
Muestra

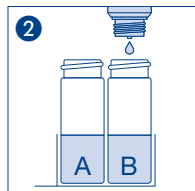


10 Δ NH_4 -1



Schütteln
Shake
Agiter
Agitar

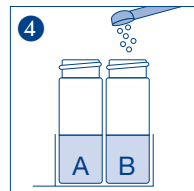


1 Δ NH_4 -2

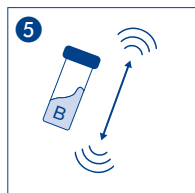


Schütteln
Shake
Agiter
Agitar




Warten
Wait
Attendre
Espere

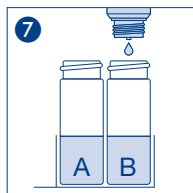


4 Δ NH_4 -3

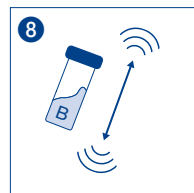

Schütteln
Shake
Agiter
Agitar

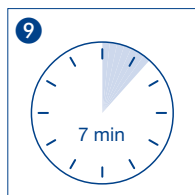

Warten
Wait
Attendre
Espere

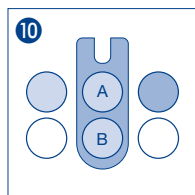

Messung
Measurement
Mesure
Medición

Geeignet für Meerwasser nach Verdünnung 1 + 9
Suitable for sea water after dilution 1 + 9
Convient à l'eau de mer après dilution 1 + 9
Adecuado para el agua de mar tras su dilución 1 + 9