



visocolor® ECO Ammonium 15

Ammonium 15 / Ammonium 15 / Amonio 15

Messbereich / Range / Domaine de mesure / Rango
0.5 – 15 mg/L NH₄⁺

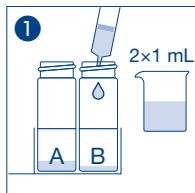


1
Probe
Sample
Echantillon
Muestra

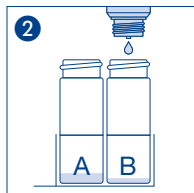


2
10 ∅ NH₄-1



3
Schütteln
Shake
Agiter
Agitar

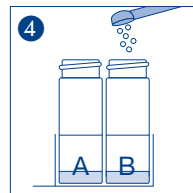


4
1 ∅ NH₄-2

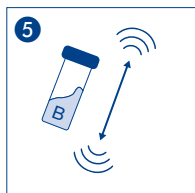


5
Schütteln
Shake
Agiter
Agitar




6
Warten
Wait
Attendre
Espere

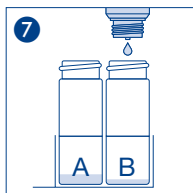


7
4 ∅ NH₄-3

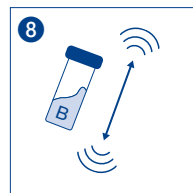

8
Schütteln
Shake
Agiter
Agitar

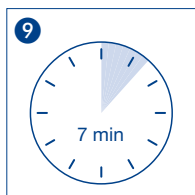

9
Warten
Wait
Attendre
Espere

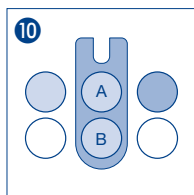

10
Messung
Measurement
Mesure
Medición

Geeignet für Meerwasser nach Verdünnung 1 + 9
Suitable for sea water after dilution 1 + 9
Convient à l'eau de mer après dilution 1 + 9
Adecuado para el agua de mar tras su dilución 1 + 9



MACHEREY-NAGEL
www.mn-net.com