



visocolor® ECO Chlorid

Chloride / Chlorure / Cloruro

Messbereich / Range / Domaine de mesure / Rango
1 – 60 mg/L Cl⁻

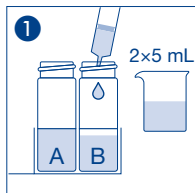


1
Probe
Sample
Echantillon
Muestra

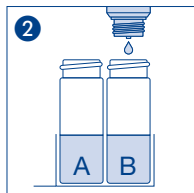


2
10 Δ Cl-1



3
Schütteln
Shake
Agiter
Agitar

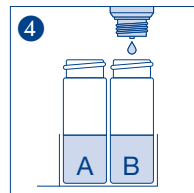


4
10 Δ Cl-2

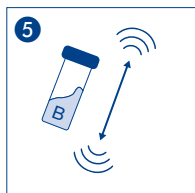

5
Schütteln
Shake
Agiter
Agitar



6
Warten
Wait
Attendre
Espere

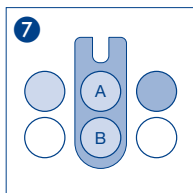

7
Messung
Measurement
Mesure
Medición

Nicht geeignet für Meerwasser
Not suitable for seawater
Ne convient pas à l'eau de mer
No es adecuado para el agua de mar